

Herculite™ *Ultra Flow*

Nanohybrid Flowable Composite



THE STRENGTH, LOW SHRINKAGE,
AND ESTHETICS YOU NEED FROM
A SINGLE FLOWABLE.



THE PERFECT DUO

For years, the Herculite™ brand has established clinical success in the practice. Now, Herculite Ultra Flow combines the long-standing expertise of the Herculite brand with a flowable composite. Plus, it's the perfect match with the esthetics and performance of Herculite Ultra.



Base/Liner
Good flow/adaptation

Case image courtesy of Dr. Bob Lowe.



Class V
Viscosity/body

Case image courtesy of Dr. Marc Geissberger.


SMART PLACEMENT TECHNOLOGY

For a flowable, having one viscosity that flows into all internal line and point angles while holding its shape may seem impossible. However, the rheological additive in Herculite Ultra Flow helps it perform this feat by reducing surface tension as the material contacts the tooth surface under pressure – allowing it to flow to the boundaries of the restoration.

The same additive encourages higher surface tension after the material has been placed. So, it won't run or drip.


A flowable with the performance you need.

OUTSTANDING FLEXURAL STRENGTH*

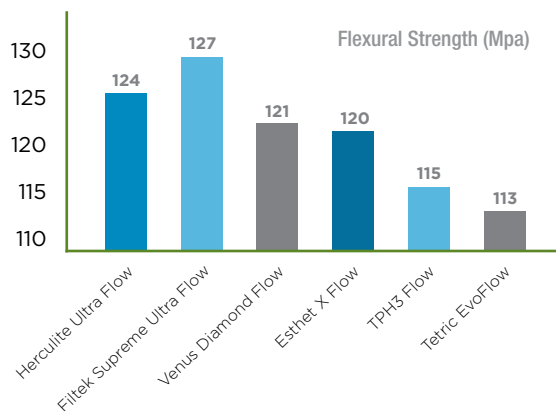


Similar to its companion Herculite™ Ultra, Herculite Ultra Flow has excellent mechanical properties, including strength.

LOW SHRINKAGE*

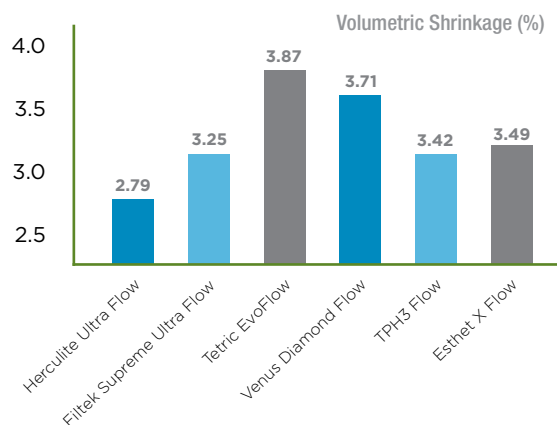


Leading experts cite low shrinkage as a critical feature of flowable composites; Herculite Ultra Flow was specially designed with a high-filler loading rate for minimal polymerization shrinkage.

OVER 300% RADIOPACITY*

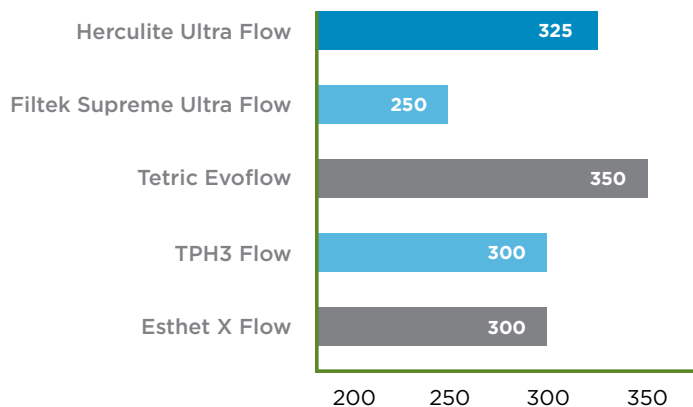

With its high radiopacity, Herculite Ultra Flow is easily discernable on radiographs.

HIGH GLOSS RETENTION*

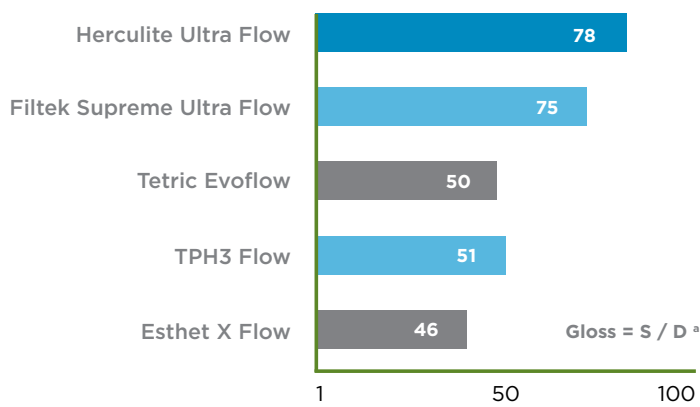

While initial gloss is a critical feature, gloss over time is more important. It is an indication of how long the esthetics of the restoration will last. Herculite Ultra Flow retains gloss after aging better than many other leading materials.

^a Contrast gloss = Specular reflection @ 60° / Diffuse reflection

ORDER TODAY!

For your convenience, Herculite™ Ultra Flow comes in 12 shades.

REFILLS:

- 2 Syringes (2g each)
- 20 Dispensing Tips
- Technique Guide
- Directions for Use

ACCESSORIES:

33370 - Tip Refills (pack of 50)



SHADES:

- 35392 A1
- 35393 A2
- 35394 A3
- 35395 A3.5
- 35396 A4
- 35397 B1
- 35398 B2
- 35399 C2
- 35400 D2
- 35401 XL1
- 35402 XL2
- 35403 Universal Opaque

RELATED PRODUCTS

BONDING AGENTS

- OptiBond™ Solo Plus™
- OptiBond™ XTR
- OptiBond™ FL
- OptiBond™ All-In-One



CURING LIGHTS

- Demi™ Ultra
- Demi™ Plus



UNIVERSAL COMPOSITE

- Herculite™ Ultra
- Premise™
- Herculite XRV™
- Point 4™



SINGLE FILL™

- SonicFill™

